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
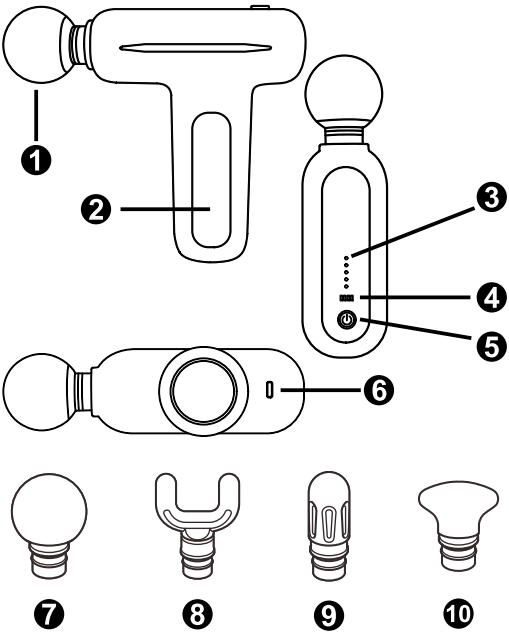
Device and controls

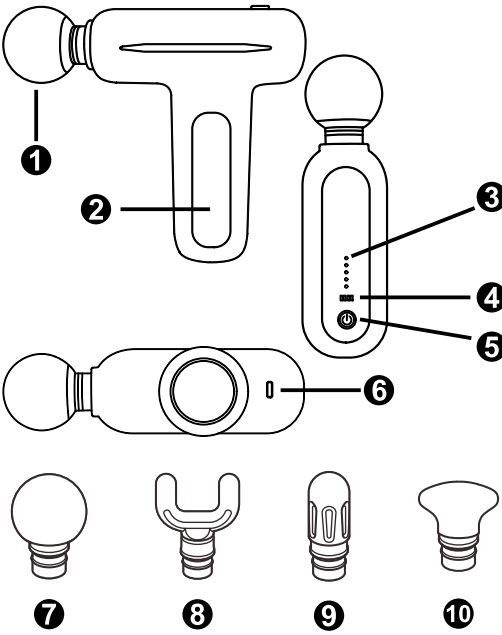Explanation of symbols

 **IMPORTANT**
Failure to follow these instructions may result in serious injury or damage to the device.

 **WARNING**
These instructions must be followed to prevent possible injury to the user.

 **CAUTION**
These instructions must be followed to prevent possible damage to the device.

 **NOTE**
These instructions provide useful additional information regarding installation and operation.

 The message device must only be used in enclosed spaces!

 Do not use the message device in damp rooms (e.g. when bathing or showering).

LOT LOT number

 Manufacturer

 Recycling symbols/codes: These are used to provide information about the material and its proper use and recycling.

 **WARNING**
Please ensure that the polythene packing is kept away from the reach of children! Risk of suffocation!

GB Safety Instructions

 **Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.**

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