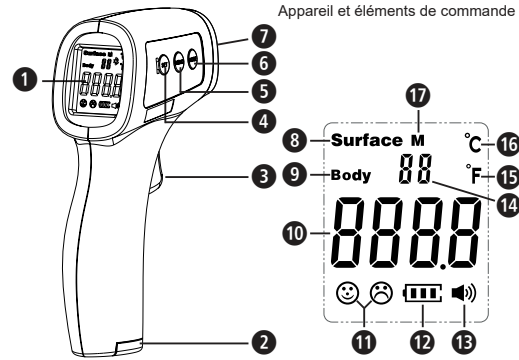
WICHTIG IMPORTANT IMPORTANT
Die Beachtung der Gebrauchsanweisung wird empfohlen.
Refer to instruction manual recommended.
Il est recommandé de respecter le mode d'emploi.

WARNING WARNING AVERTISSEMENT
Diese Warnhinweise müssen eingehalten werden, um mögliche Verletzungen des Benutzers zu verhindern. These warning notes must be observed to prevent any injury to the user. Ces avertissements doivent être respectés pour éviter d'éventuelles blessures à l'utilisateur.

ACHTUNG CAUTION ATTENTION
Diese Hinweise müssen eingehalten werden, um mögliche Beschädigungen am Gerät zu verhindern. These notes must be observed to prevent any damage to the device. Ces instructions doivent être respectées afin d'éviter d'éventuels dommages à l'appareil.

HINWEIS NOTE REMARQUE
Diese Hinweise geben Ihnen nützliche Zusatzinformationen zur Installation oder zum Betrieb. These notes give you useful additional information on the installation or operation. Ces instructions vous procurent des informations supplémentaires utiles sur l'installation ou le fonctionnement.

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Hetaida Technology Co., Ltd.

4F,5F,6F, Comprehensive Building, Baishida High-Tech Park, Ailingkan Village, Dalingshan Town, Dongguan City, Guangdong Province, P.R.China.

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GB Safety Information

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.

Intended use
The **medisana** non-contact infrared body thermometer **TM A79** is designed to be used for intermittent measurement and monitoring of human body temperature by consumers in the home setting and primary care setting as screening tool. Additionally the device is suitable to measure the temperature of objects in the range of 0°C to 100°C (32°F - 212°F).

- Only use the unit for the specific purpose described in the instruction manual. Any misuse will void the warranty.
- The unit is not designed for commercial use.
- Using this thermometer is no substitute for consulting a doctor.
- This product is designed to measure human body temperature on the forehead. Do not use it for any other purpose.
- Do not use the thermometer if it malfunctions or has been damaged in any matter.
- This product must be operated in a stable environment. If the ambient conditions are changeable, the sensor may mist up, causing an incorrect measurement. Clean the sensor according to the section „Care and Storage“.
- When the ambient temperature of the thermometer changes too much, such as moving the Thermometer from one place of lower temperature to another place of higher temperature, allow the thermometer to remain in a room for 30 minutes where the temperature is between 15°C to 40°C (59°F - 104°F).
- Remove batteries if equipment is not likely to be used for long time.
- This product is not waterproof**, do not immerse in water or other liquid; for cleaning and disinfection, please follow the „Care and Storage“ section requirements.
- Do not touch the IR sensor with your fingers.
- Avoid taking fever measurements if measures have been taken beforehand to cool the user (e.g. cold forehead compresses, etc.).
- For the measurement of human forehead temperature, select „body“ mode; for measuring other objects, liquids, food and other temperatures, select „surface“ mode.
- Do not use near strong electrostatic or magnetic fields to avoid inaccurate measurement results.
- It may affect the accuracy of measurements when the forehead is covered by hair, perspiration, cap or scarf.
- Children must be supervised to ensure that they do not play with the device.
- The temperature in the forehead and temple area differs from the internal temperature, which is taken orally or rectally. Vasoconstriction, an effect which constricts the blood vessels and cools the skin, can occur during the early stages of a fever. In this case, the temperature measured by this thermometer may be unusually low. If the measurement therefore does not match the patient's own perception or is unusually low, repeat the measurement every 15 minutes. As a reference, you can also measure the internal body temperature using a conventional oral or rectal thermometer. A control measurement using a conventional thermometer is recommended in the following cases:
 - If the reading is surprisingly low,
 - For new-born infants up to 100 days old,
 - For children under three years of age who have a weakened immune system or who react unusually in the presence or absence of fever,
 - If you are using the thermometer for the first time or want to familiarise yourself with it.
- Babies' skin reacts very quickly to the ambient temperature. Therefore, do not take their temperature with this thermometer during or after breastfeeding, because the skin temperature may then be lower than the internal body temperature.
- The device should be kept out of the reach of children/pets.
- When not in use, store the device in a dry room and protect it against extreme moisture, heat, lint, dust and direct sunlight. Never place any heavy objects on the device.
- This thermometer will not replace diagnostics in hospitals.
- Do not disassemble or modify the device. Do not let it fall down.
- Do not use this device if you think it is damaged or if you notice anything unusual. The continuous use of a damaged unit may cause injury, improper results or other serious dangers.
- This device comprises sensitive components and must be treated with caution. Observe the storage and operating conditions described in the „Technical Specifications“ section.
- No servicing/maintenance while the thermometer is in use.
- Never touch batteries and the patient simultaneously.
- Based on the current science and technology, no potential allergic reactions are known.
- This equipment needs to be installed and put into service in accordance with the information provided in this instruction manual.

Items supplied and packaging

Please check first of all that the device is complete and is not damaged in any way. In case of doubt, do not use it and contact your dealer or your service centre.

The following parts are included:

- 1 **medisana** infrared body thermometer **TM A79**
- 2 batteries 1.5V (type AAA/LR03)
- 1 Instruction manual

The packaging can be reused or recycled. Please dispose properly of any packaging material no longer required. If you notice any transport damage during unpacking, please contact your dealer without delay.

WARNING - Pay attention, that the polythene packing is kept away from the reach of children! Risk of suffocation!

What is fever?

- Body temperature varies from person to person and fluctuates during the course of the day. For this reason, it is suggested to know one's normal, healthy forehead temperature to correctly determine the temperature.
- Body temperature runs approximately from 35.5°C to 37.8°C (95.9°F-100°F). To determine if one has a fever, compare the temperature detected with the person's normal temperature. A rise over the reference body temperature of 1°C (1°F) or more is generally indication of fever.
- Different measurement sites (rectal, axillary, oral, frontal, auricular) will provide different readings. Therefore it is wrong to compare the measurement taken from different sites.
- Here below are typical temperatures for adults, based on different measurement sites:
 - Rectal: 36.6°C - 38°C / 97.9°F - 100.4°F
 - Axillary: 34.7°C - 37.3°C / 94.5°F - 99.1°F

Unit and controls

- Backlight LCD
- Battery compartment
- ON-button (Start of measurement)
- SET-button
- MEMO-button
- MODE-button
- IR Sensor
- Surface temperature measurement mode
- Body temperature measurement mode
- Measurement result
- Indicator of measurement result
- Low battery indicator
- Acoustic signal symbol
- Memory display
- Measurement unit: ° Fahrenheit
- Measurement unit: ° Celsius
- Memory indicator

Removing / inserting batteries

If the low battery indicator 12 or nothing appears in the LC-display 10 after a certain period of use, you must replace the batteries with two new 1.5V-batteries (type AAA / LR03).

Battery replacement: Slide the lid of the battery compartment cover downwards and take the cover off. Take out the batteries. Put in two new 1.5V-batteries (AAA / LR03). Make sure the batteries are the right way around. Put the battery compartment cover back on and move it upwards until it locks into position.

- Do not disassemble batteries!
- Clean the battery and device contacts if necessary before putting in the batteries!
- Remove discharged batteries from the device immediately!
- Do not use rechargeable batteries!
- Increased risk of leakage, avoid contact with skin, eyes and mucous membranes!
If battery acid comes in contact with any of these parts, rinse the affected area with copious amounts of fresh water and seek medical attention immediately!
- If a battery has been swallowed seek medical attention immediately!
- Insert the batteries correctly, observing the polarity!
- Keep the battery compartment well sealed!
- Remove the batteries from the device if it is not going to be used for an extended period!
- Keep batteries out of children's reach!
- Do not attempt to recharge the batteries! There is a danger of explosion!
- Do not short circuit! There is a danger of explosion!
- Do not throw into a fire! There is a danger of explosion!
- Keep unused batteries in their packaging away from metal objects in order to prevent short circuiting!

The ambient temperature should be stable. Keep device away from the larger flow fans, direct sunlight, air-conditioning vents and so on. Wait at least 30 seconds between each measurement. If for some reason the forehead temperature measurement cannot be realized, try to align behind the ears (see below).

When the device is running, pressing the **MODE**-button **6** to switch between „body“ mode and „surface“ mode. „Body“ mode is used for measuring human body temperature, the „Surface mode“ is used to measure the surface temperature (the factory default is „body“ mode).

- Turn on the thermometer by pressing the **ON-button** ③. The thermometer will perform a self-test with all segments displayed for 2 seconds.
- Immediately align the thermometer with the middle of the forehead to measure body temperature (between the eyebrows) and keep a distance of 1 - 5 cm. Alternatively you may use an alignment behind the ears. After measurement, a sound will be emitted and the temperature will be displayed. If you hear three beeps, the measurement was not successful. Press the **ON-button** ③ again and repeat the measurement.
- Depending on the measured body temperature the backlight LCD display will change the colors: Green means normal temperature (35.5°C-37.3°C / 95.9°F-99.1°F), yellow means slight fever (37.4°C-38.0°C / 99.3°F-100.4°F) and red means high fever (38.1°C-43.0°C / 100.6°F-109.4°F).

1. Turn on the thermometer by pressing the **ON**-button **3**. The thermometer will perform a self-test with all segments displayed for 2 seconds.
2. Aim the sensor at the center of the object at a distance of min. 5 to max. 15 cm. If it is covered with water, dust or dirt, remove this first to optimize the measurement result.
3. In surface mode the backlight is always green.

The last temperature taken before the thermometer powers off is automatically stored in memory (max. 50 results). With the unit switched off, press the **MEMO-button** to call up the values stored in memory. If not all 50 memory locations are occupied, "----C" or "----F" may appear. If you press and hold the **MEMO-button** for about 2 seconds during the memory display, "**CLR**" appears in the display and the stored data is deleted. In addition, an acoustic signal sounds.

The device allows some individual settings (F1 to F3): Press and hold the **SET**-button when the instrument is switched on to modify these settings.

F1: Measurement unit. °C or °F.
F2: Fever alarm temperature
F3: Sound on/off

With the device switched on, press and hold the **MODE**-button **6** until „rst“ appears on the display. 2 seconds later the factory settings for parameters F1 to F3 are restored (F1=°C; F2=38.1°C; F3=on).

Display	Problem	Solution
°C 	Measured body temperature outside normal limits (34.0 - 43.0°C / 93.2 - 109.4°F).	Repeat the measurement in „Body“ mode in the correct way.
°C 		
°C 	Operating temperature not correct	Establish the correct operating conditions and wait 30 minutes before operating the device.
°C 	Distance more than 1-5 cm	The optimal distance is 1 cm
	Incorrect measuring position	Please observe the instructions
	Sweat, forehead hairs or similar obstacles	Remove the obstacles
 	Display flickers or turns off	Replace batteries. The device may also be damaged. Contact the service center.
	Ambient temperature changes too quickly	Environmental conditions must be stable
	Battery capacity too low	Replace batteries.
 (blank display)	Power is off / Batteries not inserted correctly / Batteries exhausted	Replace the batteries. Insert the batteries correctly. Turn on the machine.

Care and Storage

- The IR sensor **②** is the most delicate part of the **medisana** infrared body thermometer **TM A79**. Protect it from dirt and damage. Clean the sensor with a cotton bud dipped in warm water or 70% isopropyl alcohol.
- Use a clean, soft cloth to clean the surface of the unit and the display. Do not use solvents or immerse the product in water or other liquids.
- Always keep the thermometer within the storage temperature range (- 20°C - +55°C or -4°F - +131°F) and humidity range (≤ 93%, non-condensing).
- It is recommended to store the thermometer in a dry, dust-free place. Do not expose the thermometer to direct sunlight, high temperature/humidity or extreme environments, as this may affect its function.

This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner.

Please remove the batteries before disposing of the device/unit. Do not dispose of old batteries with your household waste, but at a battery collection station at a recycling site or in a shop. Consult your municipal authority or your supplier for information about disposal.

The current version of this instruction manual can be found under www.medisana.com

The CE marking of the device refers to the EU Directive 93/42/EEC. The device is designed for domestic use. The thermometer is initially calibrated at the time of manufacture. If this thermometer is used according to the instruction manual, periodic re-adjustment is not required. If you should question the accuracy of temperature measurements, please contact the service center. Electromagnetic compatibility: see separate enclosure.

Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is. The following warranty terms apply:

1. The warranty period for medisana products is three years from date of purchase.
In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.

2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. All damage which has arisen due to improper treatment, e.g. nonobservance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

Share Info Consultant Service LLC
 Repräsentanzbüro, Heerdt Lohweg 83, 40549 Düsseldorf
 Tel: 0049 1767 0057022, Fax: 0049 1767 0057022
 E-mail: eu-rep@share-info.cn; info@share-info.cn; will.chan@share-info.cn

The service centre address is shown on the attached leaflet.

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